

Bard Orbiter™

16- & 24-pole woven diagnostic catheter

The woven advantage for simplified atrial mapping

Orbiter woven atrial mapping catheters offer exceptional stability, easy positioning and superb signal quality. Bard's unique woven material allows the Orbiter catheter to soften at body temperature to conform to the cardiac wall. The result is excellent wall contact and crisp, clear electrograms.

The Orbiter catheter provides the comprehensive coverage needed to localize the origin of most atrial tachycardias.


Exceptional stability against the atrial wall enables the Orbiter catheter to deliver consistently crisp electrograms.

A choice of curves accommodates variations in your patients' atrial anatomy.

- Woven, multi-pole construction assures reliable performance in a wide variety of clinical settings.
- Woven bioflexibility assists in positioning the curve for recording along the septum, roof and lateral wall of the right atrium.
- Bard offers the Orbiter catheter in standard and higher-torque models to fit a wide range of user preferences.



The Orbiter catheter curves cover the most common ectopic foci, allowing rapid identification of the site of earliest activation.

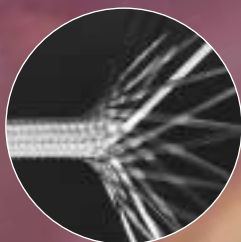

To produce woven bioflexibility, the woven core is impregnated with a proprietary polymer in an environment of strictly-controlled temperature.



Orbiter-Small, Orbiter-Large and Orbiter Milwaukee curves provide outstanding diagnostic flexibility.

Ordering Information

Bard Standard 16-Pole Orbiter™ Woven Catheters

Item No.	French	No. of Electrodes	Curve	Electrode Spacing	Tip Electrodes	Proximal Electrodes	Recording Span	Usable Length
161601	6F	16	Small 60	1,5,1mm	2mm	1mm	6.0cm	110cm
161602	6F	16	Small 60	1,5,1mm	2mm	1mm	6.0cm	80cm
161603	6F	16	Damato	1,5,1mm	2mm	1mm	6.0cm	80cm

Bard Standard 24-Pole Orbiter Woven Catheters

Item No.	French	No. of Electrodes	Curve	Electrode Spacing	Tip Electrodes	Proximal Electrodes	Recording Span	Usable Length
242401	6F	24	Milwaukee	2,4,2mm	2mm	1mm	9.3cm	110cm
242402	6F	24	Orbiter-Small	2,4,2mm	2mm	1mm	9.3cm	110cm
242403	6F	24	Orbiter-Small	2,7,2mm	2mm	1mm	12.6cm	110cm
242404	6F	24	Orbiter-Large	2,7,2mm	2mm	1mm	12.6cm	110cm
242405	6F	24	Orbiter-Small	2,9,2mm	2mm	2mm	17.1cm	110cm
242406	6F	24	Orbiter-Large	2,9,2mm	2mm	2mm	17.1cm	110cm

Bard 24-Pole High-Torque Orbiter Woven Catheters

Item No.	French	No. of Electrodes	Curve	Electrode Spacing	Tip Electrodes	Proximal Electrodes	Usable Length
242411	6F	24	Orbiter-Small	2,4,2mm	2mm	1mm	110cm
242412	6F	24	Orbiter-Small	2,7,2mm	2mm	1mm	110cm
242413	6F	24	Orbiter-Large	2,7,2mm	2mm	1mm	110cm
242414	6F	24	Orbiter-Small	2,9,2mm	2mm	1mm	110cm
242415	6F	24	Orbiter-Large	2,9,2mm	2mm	1mm	110cm
242416	6F	24	Milwaukee	2,4,2mm	2mm	1mm	110cm

Bard Orbiter Cable Information

Item No.	Use	Length	Description
200774P	Single Use	125cm	Protected 8-Pin Easy-Mate™

Bard Electrophysiology Division
USA
Tel: 978-441-6202
800-824-8724 (toll free w/in the US)
Fax: 978-323-2222
877-281-1492 (toll free w/in the US)

Bard Asia Pacific Division
USA
Tel: 908-277-8000
Fax: 908-277-8090

Bard Australia, Pty Ltd
AUSTRALIA
Tel: 61-2-8875-4000
Fax: 61-2-8875-4050

Bard Benelux, N.V.
BELGIUM
Tel: 32-0-1428-6950
Fax: 32-0-1428-6955

Bard Canada, Inc
CANADA
Tel: 905-275-8000
Fax: 905-275-3383

Bard International, Inc
CHINA
Tel: 86-10-8391-3282
Fax: 86-10-8391-3285

Laboratories Bard SA
FRANCE
Tel: 33-1-39-30-5858
Fax: 33-1-39-30-5822

C. R. Bard Germany, GmbH
GERMANY
Tel: 49-721-9445-265
Fax: 49-721-9445-111

Bard International
INDIA
Tel: 91-22-612-3783
Fax: 91-22-612-3821

Bard S.p.A.
ITALY
Tel: 39-06-52-493-1
Fax: 39-06-52-958-52

Medicon, Inc
JAPAN
Tel: 6-6203-6546
Fax: 6-6203-1516

Bard International, Inc.
KOREA
Tel: 82-2-719-3080
Fax: 82-2-719-3081

Bard Latin America
USA
Tel: 908-277-8000
Fax: 908-277-8090

Bard Medical Systems Norden
SWEDEN
Tel: 46-42-386-000
Fax: 46-42-386-010

C. R. Bard Portugal, LDA
PORTUGAL
Tel: 351 21 319 03 30
Fax: 351 21 319 03 39

Bard Singapore Pte. Ltd.
SINGAPORE
Tel: 65-250-7311
Fax: 65-353-7311

Bard De Espana, SA
SPAIN
Tel: 34-93-253-7800
Fax: 34-93-253-7834

Bard Medica, SA
SWITZERLAND
Tel: 41-22-885-9618
Fax: 41-22-885-9635

Pacific Health Care Co., Ltd
TAIWAN
Tel: 886-2274-07-891
Fax: 886-2277-75-576

Bard Limited
UNITED KINGDOM
Tel: 44-1-293-527-888
Fax: 44-1-293-552-428



Bard Electrophysiology
55 Technology Drive • Lowell, MA 01851 USA
Tel: (978) 441-6202 • Fax: (978) 323-2222
www.BARDEP.com



Before use, consult product labels and inserts for any indications, contraindications, hazards, warnings, cautions and use.

Bard and the stylized heart are registered trademarks of C.R. Bard, Inc. or an affiliate. Orbiter and Easy-Mate are trademarks of C.R. Bard, Inc. or an affiliate. All trademarks or registered trademarks are the property of their respective owners.
© 2001 C.R. Bard, Inc. All rights reserved. LT04Z0004/Rev. 00/06-01