

# Orbiter® PV

Variable Loop Mapping Catheter

We've Got  
Contact



# Orbiter® PV

## Variable Loop Mapping Catheter

# We've

When it comes to pulmonary vein isolation in the treatment of focal atrial fibrillation, the challenges can be numerous. The Orbiter PV Variable Loop Mapping Catheter has been designed to overcome these challenges. With a reduction of complications, fluoroscopy and procedural time in mind, we've created a catheter that keeps you in contact from mapping through ablation and confirmation of post-ablation conduction block.

### Location is Everything

It's no secret that true ostial positioning is a key element in successful and lasting outcomes. The Orbiter PV catheter is the first circular mapping catheter to combine an **adjustable** loop with a **deflectable** shaft. The variable loop design allows you to adjust the catheter's diameter to maximize contact with the pulmonary vein at its most proximal point.

### Optimal Signal Quality

The clear signal quality that you are accustomed to from Bard Electrophysiology is only enhanced by the ability to adjust loop size and maximize electrode-to-tissue contact with outward radial pressure.



As the loop expands, radial pressure is extended outward providing stability and securing the position



Variable Loop  
Thumbwheel

Deflectable Shaft Slider

# Got Contact


## Stability

One of the greatest challenges faced by electrophysiologists today is establishing ostial position and maintaining orientation throughout the procedure. The Orbiter PV product is built around a pressure-based positioning feature that maximizes stability in the physician-preferred location.

## Procedural Efficiency

Among individual patients, pulmonary vein anatomy varies widely. By incorporating the ability to adjust loop diameter from 14.5mm to 25mm, we reduce the need for multiple fixed-diameter mapping catheters and their associated exchanges.

For more information about The Orbiter PV product, call your Bard Representative or visit our website at [www.bardep.com](http://www.bardep.com)

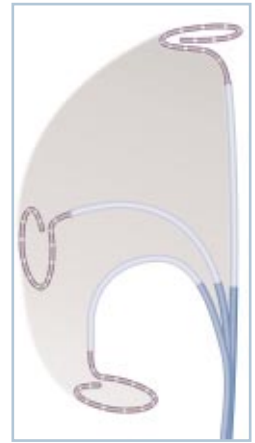

The deflectable shaft offers a complete range of motion



With continuous thumbwheel-controlled precision, the Orbiter PV catheter expands and contracts to fit virtually any anatomical situation



# Ordering Information

## Orbiter PV Variable Loop Mapping Catheter

Catalog Number	French Size	Expandable Loop	No. of Electrodes	Electrode Spacing	Usable Length
320100	7.5F (5F tip)	14.5mm to 25mm	14	5mm	105mm

Use with 8F Introducer

## Cable Information

Item No.	Description	Length
390003P	Protected 24-pin cable	126cm
390004P	Protected 24-pin cable	210cm

### Bard Electrophysiology Division

USA  
Tel 978-441-6202  
800-824-8724 (toll free w/in the US)  
Fax 978-323-2222  
877-281-1492 (toll free w/in the US)

### Bard Asia Pacific Division

USA  
Tel 908-277-8000  
Fax 908-277-8090

### Bard Australia, Pty Ltd

AUSTRALIA  
Tel 61-2-8875-4000  
Fax 61-2-8875-4050

### Bard Medica S.A.

AUSTRIA  
Tel 43-1-4949-130  
Fax 43-1-4949-13030

### Bard Benelux, N.V.

BELGIUM  
Tel 32-0-1428-5900  
Fax 32-0-1428-6955

### Bard Canada, Inc.

CANADA  
Tel 905-275-8000  
Fax 905-275-5370

### Bard International, Inc.

CHINA  
Tel 86-10-6501-4225  
Fax 86-10-6501-2573

### Laboratories Bard SA

FRANCE  
Tel 33-1-39-30-5858  
Fax 33-1-39-30-5859

### C. R. Bard Germany, GmbH

GERMANY  
Tel 49-721-9445-265  
Fax 49-721-9445-111

### Bard South Asia

INDIA  
Tel 91-22-612-3783  
Fax 91-22-612-3821

### Bard S.p.A.

ITALY  
Tel 39-06-52-493-1  
Fax 39-06-52-950-67

### Medicon, Inc.

JAPAN  
Tel 6-6203-6546  
Fax 6-6203-1516

### Bard International, Inc.

KOREA  
Tel 82-2-719-3080  
Fax 82-2-719-3081

### Bard Latin America

USA  
Tel 908-277-8000  
Fax 908-277-8090

### C. R. Bard Portugal, LDA

PORTUGAL  
Tel 351-21-319-0330  
Fax 351-21-319-0339

### Bard International, Inc.

SINGAPORE  
Tel 65-250-7311  
Fax 65-353-7311

### Bard De Espana, SA

SPAIN  
Tel 34-93-253-7800  
Fax 34-93-205-8229

### Bard Medical Systems Norden

SWEDEN  
Tel 46-42-386-000  
Fax 46-42-386-010

### Bard Medica, SA

SWITZERLAND  
Tel 41-17-225-360  
Fax 41-17-225-370

### Pacific Health Care Co., Ltd

TAIWAN  
Tel 886-2274-07-891  
Fax 886-2277-75-576

### Bard Limited

UNITED KINGDOM  
Tel 44-1-293-527-888  
Fax 44-1-293-552-428

### Bard Sdn Bhd

WEST MALAYSIA  
Tel 60-4-489-1300  
Fax 60-4-489-1737



Bard Electrophysiology Division  
C. R. Bard, Inc.  
55 Technology Drive • Lowell, MA 01851 USA  
Tel: (978) 441-6202 • Fax: (978) 323-2222  
[www.BARDEP.com](http://www.BARDEP.com)



Before use, consult product labels and inserts for any indications, contraindications, hazards, warnings, cautions and use.

Bard and the stylized heart and Orbiter are registered trademarks of C.R. Bard, Inc. or an affiliate. All trademarks or registered trademarks are the property of their respective owners.  
© 2002 C.R. Bard, Inc. All rights reserved. LT04Z0035/Rev. 01/03-03